



DR50-A01



DR90-A01

極致的電動代步車駕駛體驗


最新的 R 系列控制器打造更勝以往的駕駛體驗，展現創新完美的操控性能，並提供動力需求與空間限制的絕佳平衡。

- 50 及 90 A 的機型隨時提供您所需的動力。
- 多不勝數的創新設計包含兩組行駛檔案、煞車燈及倒車燈、倒退喇叭與電磁煞車釋放開關。
- 可調整的加速曲線、坡道上零倒滑、精密的馬達匹配演算程式與 SRW 功能，提供了最佳的障礙跨越與上坡起步的能力，帶給您極致的駕駛體驗。
- 智慧型的馬達以及電池管理特色提供最佳化的動力驅動、可選擇多種預設的電池管理程式、5 伏特及 12 伏特電量表的輸出並且深入地紀錄與分析電池的使用情形。
- 先進的診斷及售後服務工具，包含了行駛紀錄及貼心的定期保養叮嚀。
- 使用工業標準的插頭與簡單的安裝方式，提供立即的系統更換。支援多種電池種類、多功能的邏輯腳位與極富彈性的 inhibit 設定。


R-series



DR50-A01

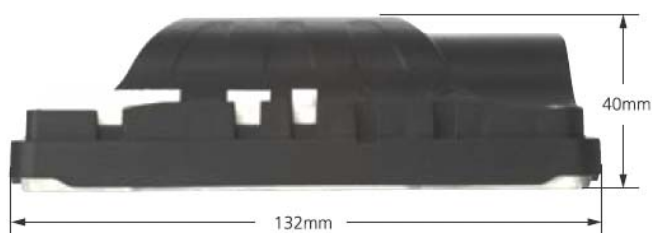

DR90-A01

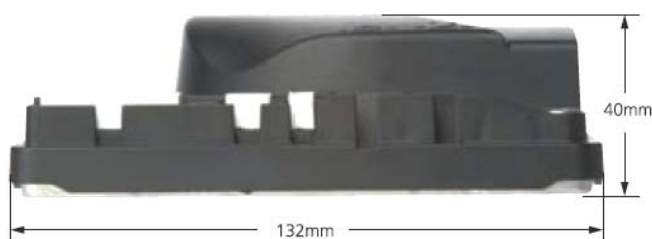

供給電壓		24 伏特 (直流)
操作電壓		18-32伏特 (直流)
最大電流	R50	40 安培
	R90	70 安培
加強電流	R50	10 安培
	R90	20 安培
瞬充時間		1-10 秒
煞車器電壓		24 伏特 (直流)
充電電流		8 安培RMS

保護等級		電子類 IPX5
操作溫度		-25°C to 50°C
保存溫度		-40°C to 65°C
材質		壓鑄鋁底部，塑膠上蓋
型號	R50	DR50-A01
	R90	DR90-A01

New Zealand (Head Office)
 Ph: +64-3-962-2519
 Fax: +64-3-962-2966
 Email: info@dynamiccontrols.com

Asia
 Ph: +886-955-335-243
 Fax: +886-225-981-562
 Email: asia@dynamiccontrols.com

USA
 Ph: +1-440-979-0657
 Fax: +1-440-979-1028
 Email: usa@dynamiccontrols.com

Europe
 Ph: +44-1905-772-321
 Fax: +44-1905-827-520
 Email: europe@dynamiccontrols.com



Dynamic Controls is the world's leading manufacturer of electronic controls for power wheelchairs and scooters. Dynamic was established in 1972 and is headquartered in New Zealand. Regional centres are located in Europe, United States, Asia and Australasia. ISO 13485 Accredited – Dynamic goes above and beyond industry standard expectations to ensure customers receive the best products possible.

For more information about Dynamic Controls visit www.dynamiccontrols.com